

## **Vertical Saturated Steam Sterilizers**




- Vertical autoclaves **DLVS** and **DLVD series** feature a simple construction and a user friendly interface. They are the ideal solution for all laboratories, pilot plants and research centres applications.
- **DLVS series** (single wall) are mainly used for solids that can be wet/ not completely dry after sterilisation and for liquid solutions in sealed/ unsealed glass containers
- **DLVD series** (with jacket and liquid ring vacuum pump) are mainly used for solids that must be dry after sterilisation and for porous materials

- Main AISI 316L stainless steel construction
- Self steam generation
- Cylindrical wire baskets for loading / unloading
- Design, manufacturing and test according to European Pressure Equipment Directive (P.E.D.) No.97/23/EC



## **Control System**

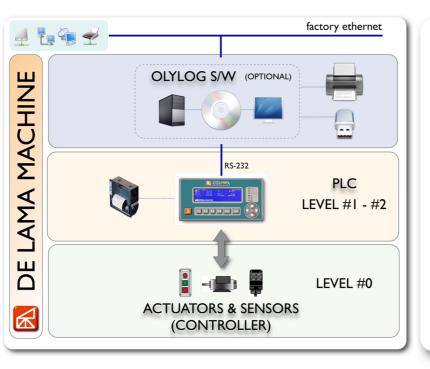

- MI and integrated PLC system WinPlus model
- Integrated UPS for power back-up
- Integrated industrial solid state memory for operating system and data storage
- Easy networking capability through optional PC based datalogger (S/W OLYLOG)

## **Available Sizes**

MODEL	CHAMBER VOLUME	CHAMBER INTERNAL DIMENSIONS Diameter Height	WEIGHT	INSTALLED POWER
	Its.	mm. mm.	kg.	kW
DLVS 40.50	71	400 500	85	4,2
<b>DLVS 40.60</b>	83	400 600	100	4,2
DLVS 40.70	95	400 700	110	4,2
DLVS 50.60	132	500 600	145	7,5
DLVS 50.70	151	500 700	165	7,5
DLVS 50.80	170	500 800	180	7,5

MODEL	CHAMBER VOLUME	INTERNAL DIMENSIONS Diameter Height		WEIGHT	INSTALLED POWER
DLVD 40.50 DLVD 40.60 DLVD 40.70 DLVD 50.60 DLVD 50.70	1ts. 71 83 95 132 151	mm. 400 400 400 500 500	mm. 500 600 700 600 700	kg. 135 155 170 215 245	5,4 5,4 5,4 9,0 9,0
DLVD 50.80	170	500	800	270	9,0

CHAMPED







